



Experiment 2

Experiment : To find the temporary and permanent hardness of water sample by complexometric titration using standard EDTA solution.


Apparatus: Pipette, burette, beakers, conical flask, funnel, burette stand and clamp.

Chemicals: water sample, ethylenediaminetetraacetic acid (EDTA), Eriochrome black-T (EBT) indicator, ammonium hydroxide - ammonium chloride buffer of pH 10.


Chemical structures

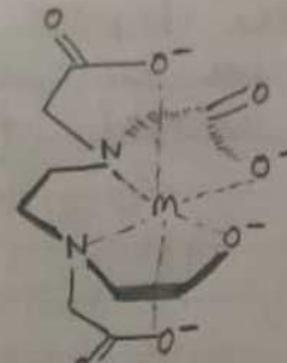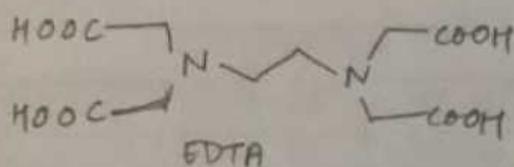

Metal EDTA complex.

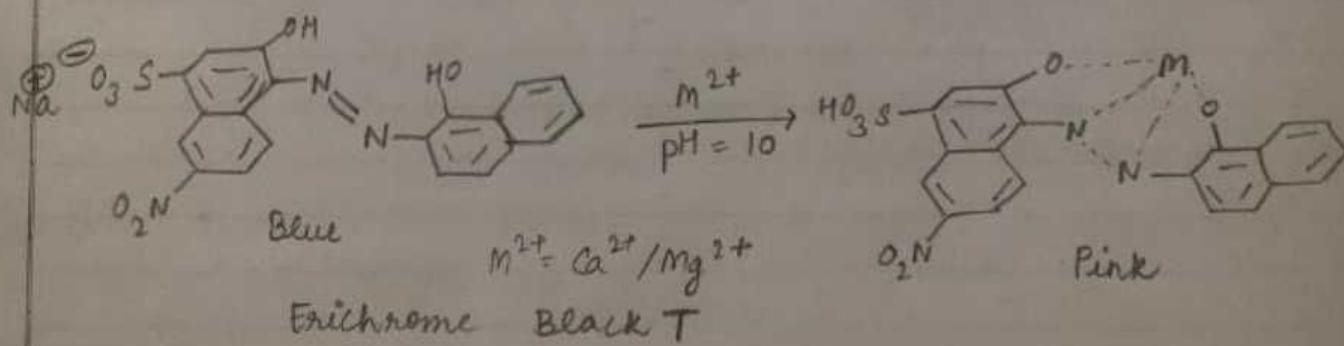

Chemical Equations

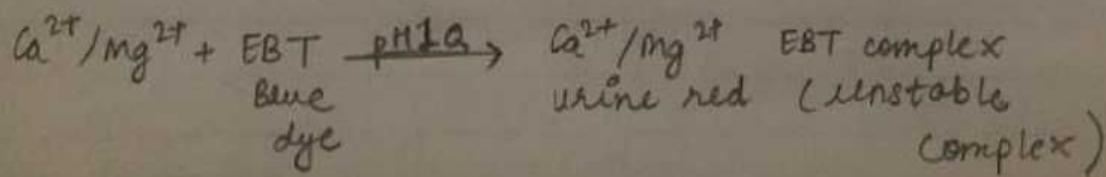


Experiment:- To find the temporary and permanent hardness of water sample by complexometric titration using standard EDTA solution.


Apparatus :- Pipette, burette, beaker, conical flask, funnel, burette stand and clamp.

Chemicals - Water samples, ethylenediaminetetraacetic (EDTA) Erichrome Black-T (EBT) indicator, ammonium hydroxide - ammonium chloride buffer of pH 10.

Theory : Hardness of water is due to the presence of solute salts of Ca and Mg. It is an important parameter to judge the quality of water. Determination of hardness of water by EDTA titration is a very accurate method based on the fact that when Erichrome black-T (EBT, blue dye) is added to the hard water (at about pH 10), it gives a ^{were} red coloured unstable complex with $\text{Ca}^{2+}/\text{Mg}^{2+}$ ion.

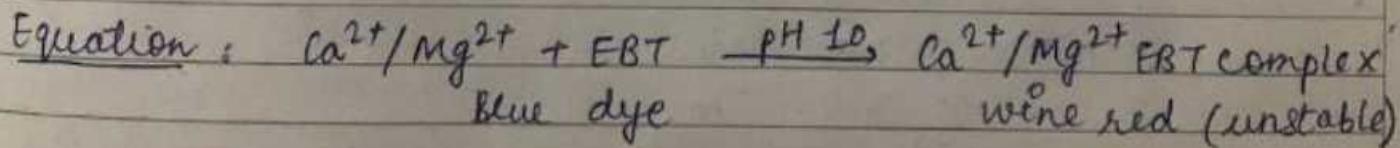

Temporary hardness in water sample is caused by bicarbonates of hardness producing ions (Ca^{2+} and Mg^{2+}). This can be removed by prolonged boiling due to decomposition of bicarbonates with the evolution of CO_2 and simultaneous precipitation of the respective carbonates, when EDTA solution is added to the hard water (with permanent and temporary

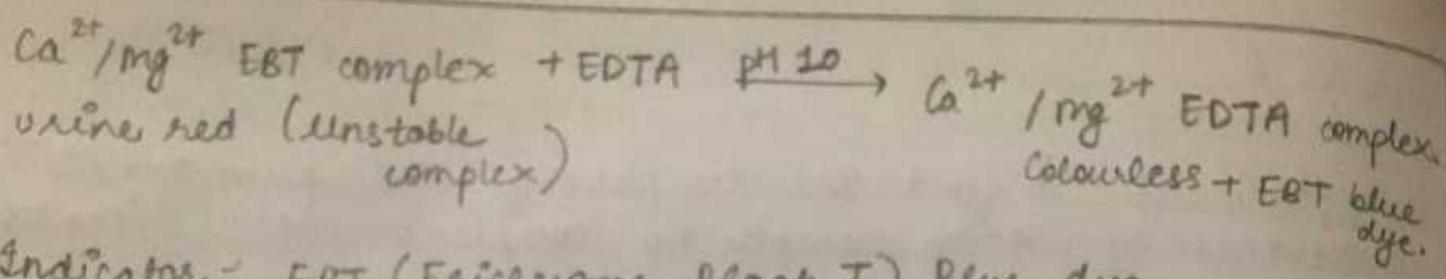

Indicator - EBT (Eriochrome Black T) Blue dye.

End point: wine red to blue colour

OBSERVATION : i) Standardisation of EDTA soln.

Volume of 0.01 M standard hard water solution taken
for each titration = 10 ml

Sr. No.	Burette reading (ml)		Vol. of EDTA used (ml)
	Initial	Final	
1	0	9.9	9.9
2	0	9.9	9.9
3	0	9.9	9.9
4	0	9.9	9.9

Mean volume of EDTA used (V_0) = 9.9 ml.

ii) Determination of total hardness


Volume of hard water sample (unknown) taken for each titration = 10ml.

S.R. No.	Burette Reading (ml)		Volume of EDTA used (ml)
	Initial	Final	
1	0	7.7	7.7
2	0	7.7	7.7
3	0	7.7	7.7
4	0	7.7	7.7

$$\text{Mean volume of EDTA used } (V_i) = \frac{7.7}{1} \text{ (ml)}$$

hardness), the unstable wine red complex of $\text{Ca}^{2+}/\text{Mg}^{2+}$ EBT breaks and a stable complex of $\text{Ca}^{2+}/\text{Mg}^{2+}$ with EDTA is formed resulting in change of color of the solution from wine red to blue at the end point.

Equation :

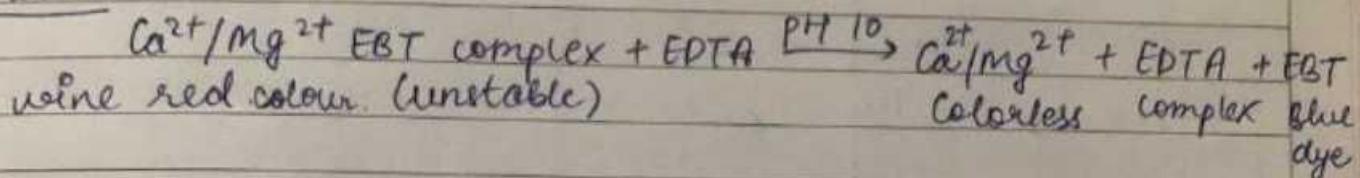

PROCEDURE :- i) Preparation of Standard Hard water:

Dissolve 1 gm of pure dry CaCO_3 in minimum quantity of dilute HCl. Evaporate the solution to dryness on a water bath to remove excess of acid. Dilute the contents with distilled water to make 1L. Each ml of the solution contains 1 mg of CaCO_3 , i.e. hardness of this solution is 1000 ppm (0.01M). This solution is used to standardize the EDTA solution.

ii) Standardization of EDTA

- a.) Rinse the titration flask with distilled water and transfer 10ml of the standard hard water sample (0.01M) into it using a pipette.
- b.) Add about 2-3 ml of ammonia/ammonium chloride buffer solution and 2-3 drops of the EBT indicator. The color of solution becomes wine red.
- c.) Titrate the hard water against the EDTA solution, till the wine red color changes to blue. Note readings (V₀ ml)
- d.) Repeat the procedure until three concordant reading are obtained.

Determination of total hardness of water sample.
Rinse the titration flask with distilled water and

iii) Determination of permanent hardness

Volume of boiled hard water sample taken for each titration = 10 ml

Sr. No.	Burette reading (ml) Initial	Burette reading (ml) Final	Vol. of EDTA used (ml)
1	0	6.3	6.3
2	0	6.3	6.3
3	0	6.3	6.3
4	0	6.3	6.3

$$\text{Mean Volume of EDTA used } (V_2) = \underline{6.3} \text{ ml.}$$

Calculations :- i) Determination of molarity of EDTA solution applying the molarity eq.

$$(\text{Standard Hard water}) = (\text{EDTA}) \\ 0.01 \times 10 = M_1 \times V_0$$

$$\text{Molarity of the EDTA } M_1 = \frac{0.01 \times 10}{V_0} \rightarrow \frac{0.01 \times 10}{9.9}$$

$$M_1 = \underline{0.0101 \text{ M}}$$

ii) Re-determination of total hardness.

Molarity of EDTA (M_1) =
Applying the molarity equation
(Hard water) = (EDTA)

$$\frac{M_2^2 \times 10}{M_2} = \frac{M_1 V_1}{10} = \frac{M_1 \times V_1}{10} \\ \frac{0.0101 \times 7.7}{10} \rightarrow 0.007777 \text{ M}$$

$$\text{Hardness of water sample, } Y = \frac{M_1 V_1}{10} \times 100 \text{ (mol wt of CaCO}_3\text{)} \text{ gm/l}$$

$$\text{Total hardness (Y)} = \frac{\frac{M_1 V_1}{10} \times 100 \times 1000 \text{ mg/L}}{777.7} \text{ mg/L} \\ \underline{777.7} \text{ ppm (mg/L)}$$

transfer 10ml of the give sample into this using pipette
 (Follow the steps 2-4 given above for standardisation of EDTA) let the titre value corresponding to total hardness of water sample be V_1 ,

Determination of permanent hardness

Make 100ml of hard water sample into 500ml beakers, boil gently for 30-35 min. Filter the solution into a 100ml measuring flask. Make up the solution upto the mark with de-ionized water and mix thoroughly. Rinse the titration flask with distilled water and transfer 10ml of this (boiled water) sample into it using a pipette. follow the step 2-4 as given above (for standardisation of EDTA) let us titre values corresponding to total hardness of water sample be V_2 .

Determination of temporary hardness

Difference between the two values ($V_1 - V_2$) corresponds to temporary hardness.

General calculations

i) Determining the molarity of EDTA soln

Applying the molarity equation

(Standard Hard water) (EDTA)

$$0.01 \times 10 = M_1 \times V_0$$

$$M_1 = (0.01 \times 10) / V_0$$

ii) Determination of total hardness

Molarity of EDTA = M_1

Applying molarity equation (Hard) (EDTA)

$$M_2 \times 10 = M_1 \times V_1$$

Molarity of hard water, $M_2 = (M_1 \times V_1) / 10$

Hardness of water sample, $y = \text{molarity} \times 100$

Teacher's Signature : _____

iii) Determination of permanent hardness.

Apply the molarity equation

(Boiled Hard water)

(EDTA)

$$M_3 \times V_1 = M_1 \times V_2 \quad (\text{EDTA})$$
$$M_3 = (\text{Molarity of hard water due to permanent hardness})$$
$$M_3 = \frac{M_1 V_2}{10} = \frac{0.0101 \times 6.3}{10} = 0.006363 \text{ M}$$

Permanent hardness of water sample (z) :

$$\frac{M_1 V_2}{10} \times 100 (\text{mol wt of } \text{CaCO}_3) \text{ gm/L}$$
$$= \frac{M_1 V_2 \times 100 \times 1000 \text{ mg/L}}{10} = \underline{\underline{636.3}}$$

$$\text{Total hardness (y)} = \underline{\underline{636.3}} \text{ ppm (mg/L)}$$

$$\text{iv) Temporary hardness} = (y-z) \text{ ppm} \Rightarrow 141.4 \text{ ppm}$$

Results \rightarrow Total hardness (y) = $\underline{\underline{777.7}}$ ppm

$$\text{Permanent hardness (z)} = \underline{\underline{636.3}} \text{ ppm}$$

$$\text{Temporary hardness (y-z)} = \underline{\underline{141.4}} \text{ ppm}$$

molecular weight of CaCO_3

Hardness of water sample, $Y = \frac{M_1 V_1}{10} \times 100$ (m.w of CaCO_3) gm/L
 $= \frac{(M_1 V_1)}{10} \times 100$ (M.W of CaCO_3) $\times 1000$ mg/L

Total hardness = Y ppm (mg/L)

Determination of permanent hardness:

(Hard Water) EDTA

$$M_3 \times V_1 = M_1 \times V_2 \quad (M_3 = \text{Molarity due to permanent hardness})$$

$$\text{Molarity of hard water} = \frac{M_1 \times V_2}{V_1} = \frac{M_1 \times V_2}{10}$$

Permanent hardness of water sample, $Z = \text{Molarity} \times \text{M.W of } \text{CaCO}_3$

Permanent hardness of water sample

$$Z = \frac{(M_1 V_2)}{10} \times 100 \text{ (Mol wt of } \text{CaCO}_3 \text{)} \text{ gm/L}$$

$$= (M_1 V_2) / 10 \times 100 \times 1000 \text{ mg/L}$$

Permanent hardness = Z ppm (mg/L)

(iv) Temporary Hardness: Total - Permanent = $(Y-Z)$ ppm.

Results: Total Hardness = Y ppm = 777.7 ppm

Permanent = Z ppm = 636.3 ppm

Temporary hardness = $(Y-Z)$ ppm = 141.4 ppm.

PRECAUTIONS \rightarrow i.) wash Titration flask with distilled water each time, before transferring hard/sample water soln.

Expected CLO & Daily life Application \rightarrow Determination of hardness of water can help in industrial settings, where water hardness is monitored to avoid costly breakdown in boilers, cooling towers and other equipment. High calcium levels also cause irritation in eyes. Hence, qualification of ions in potable water is necessary.