


Experiment 1

Experiment :- To determine the amount of NaOH and Na_2CO_3 present in the same solution.


Apparatus :- Pipette, burette, beakers, conical flask, funnel, burette stand and clamp.

Chemicals \rightarrow Sodium carbonate (Na_2CO_3), Sodium Hydroxide (NaOH), hydrochloric acid (HCl), methyl orange and phenolphthalein.

Chemical structures

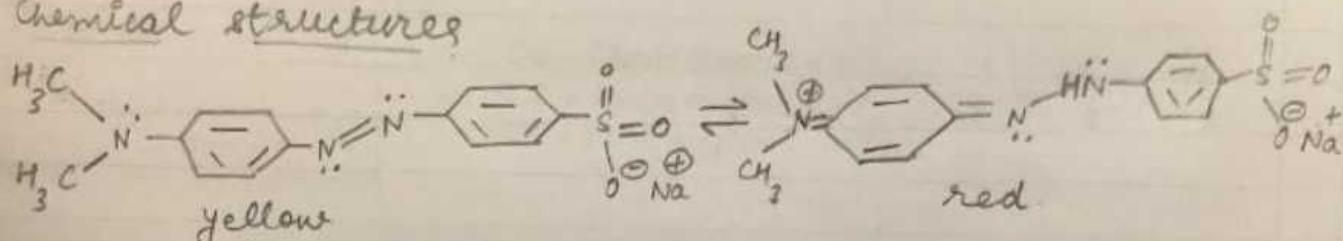


Methyl orange.

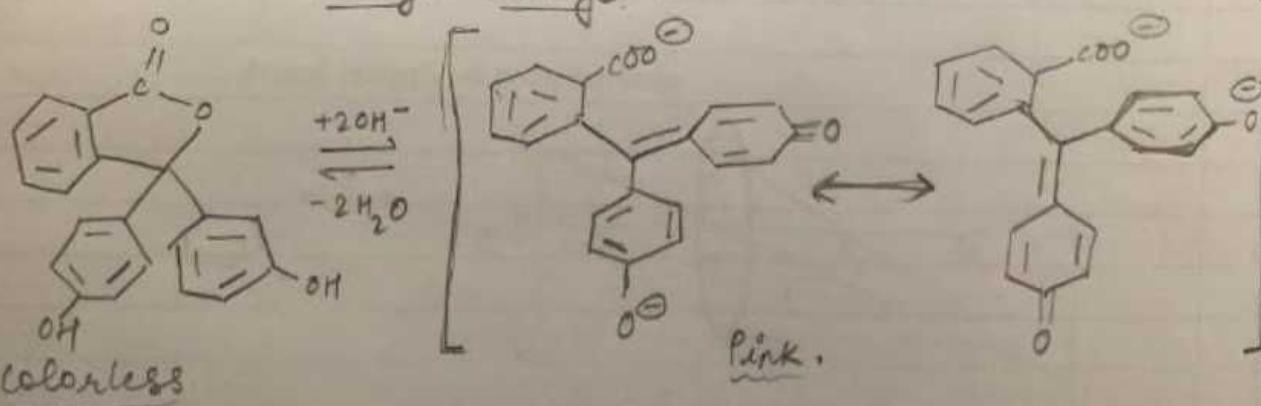

Phenolphthalein

Chemical Equation

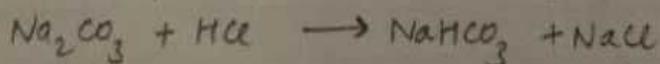


Kushalji

Experiment :- To determine the amount of NaOH and Na_2CO_3 present in the same solution.

Apparatus :- Pipette, burette, beakers, conical flask, funnel, burette stand and clamp.

Chemicals - Sodium carbonate (Na_2CO_3), sodium hydroxide (NaOH), hydrochloric acid (HCl), methyl orange and phenolphthalein

Theory :- When a mixture of NaOH and Na_2CO_3 is titrated against a standard HCl solution, color of the solution changes from yellow to pink (using methyl orange as an indicator) due to complete neutralisation of both the alkalis at $\text{pH} \approx 4$, when the mixture is titrated with phenolphthalein as an indicator, the colour of solution changes from pink to colorless due to complete neutralisation of NaOH and half neutralisation of Na_2CO_3 (i.e upto the conversion of Na_2CO_3 to NaHCO_3) at $\text{pH} \approx 8$. The difference of two titre values gives the amount of HCl required for half neutralisation of Na_2CO_3 while the difference of first value and twice the second titre values gives the amount of HCl required for NaOH neutralisation.

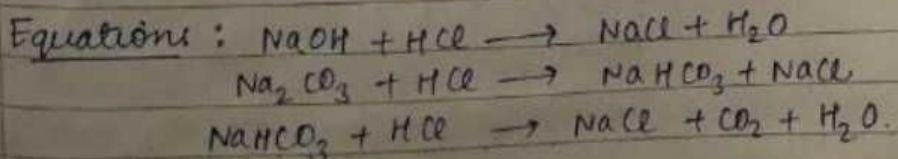

Indicator \rightarrow Methyl orange, Phenolphthalein.

end point \rightarrow Yellow to orange.

Observation: (i) standardization of HCl soln.

Volume of 0.1 N Na_2CO_3 solution taken for each titration
= 10 ml.

S. No	Burette Reading (ml)	Volume of the HCl used (ml)
	Initial	Final
1	0	9.9
2	0	9.9
3	0	9.9

$$\text{Mean volume of HCl} (V_0) = 9.9 \text{ ml}$$

(ii) Determination of NaOH and Na_2CO_3 in the mixture

Volume of mixture of NaOH and Na_2CO_3 solution taken for each titration = 10 ml.

S. No	Burette reading (ml)			Vol. of HCl used (ml)	
	Initial (A)	colorless with phenolphthalein (B)	orange color with methyl orange (C)	P = B - A	M = C - A
1	0	14 - 0.1 = 13.9	16.5 - 0.1 = 16.4	13.9	16.4
2	0	13.9	16.4	13.9	16.4
3	0	13.9	16.4	13.9	16.4

$$\text{Mean value of HCl used for P} = 13.9 \text{ ml}$$

$$\text{Mean value of HCl used for M} = 16.4 \text{ ml}$$

As P corresponds to $\frac{1}{2}$ neutralisation Na_2CO_3 and NaOH

$$\text{Half of } \text{Na}_2\text{CO}_3 = M - P = 2.5 \text{ ml}$$

X. 1000

Procedure: (i) Standardisation of HCl.

- i) Transfer 10 ml's of standard 0.1 N Na_2CO_3 solution in a clean conical flask using a pipette.
 - ii) Add 2 drops of methyl orange indicator. Titrate the solution against HCl from the burette.
 - iii) The color of the solution changes from yellow to pink (end point)
 - iv) Note the volume of the solution used and repeat the titration at least 4 times and take the mean of the closely related readings (V_a)
- (ii) Determination of NaOH and Na_2CO_3 present.
- i) Transfer 10 ml of mixture of alkali solution into a conical flask.
 - ii) Add 2-3 drops of phenolphthalein indicator. The solution becomes pink in color.
 - iii) Note the initial reading of HCl from the burette (A). Titrate the solution with standard HCl while the solution becomes colourless (B).
 - iv) Note the titre values and this is the phenolphthalein end point (P). To the same solution, add 2-3 drops of methyl orange indicator and continue the titration with HCl, initial a sharp color change occurs from yellow to orange at the end point (C).
 - v) This titre value i.e., the total value to HCl run down from the beginning of the experiment to the methyl orange end point is noted and this is the methyl orange end point (M).

So, volume of HCl required for neutralization of Na_2CO_3

$$\Rightarrow 2(M-P) = V_1 = 2 \times 2.5 = 5 \text{ ml.}$$

Volume required for neutralisation of NaOH

$$V_2 = M - 2(M-P)$$

$$\Rightarrow 16.4 - 5 \Rightarrow 11.4 \text{ ml}$$

Calculation: (i) Standardization of HCl.

Vol. of alkali solution (Na_2CO_3) taken = 10 ml.

Normality (N_1) of $\text{Na}_2\text{CO}_3 = 0.1$; Vol. of HCl used = 9.9 ml (V_1)

$$N_1 = \frac{N_2}{V_0} \times 10 = \frac{0.1}{9.9} \times 10 \Rightarrow \frac{1}{9.9} = 0.101 \text{ N.}$$

(ii) Determination of NaOH and Na_2CO_3 .

(iii) Determination of NaOH

1L of 1N HCl = 40g of NaOH

Normality of HCl used = $N_1 = 0.101 \text{ N}$

$$11.4 \text{ ml of } 0.101 \text{ N HCl} = 40 \times \frac{11.4 \times 0.101}{1000} = y_1 \text{ gram of NaOH}$$
$$= 0.046 \text{ g.}$$

$$\text{Strength of NaOH} = y_1 \times 100 \Rightarrow 0.046 \times 100 = 4.6 \text{ g/L}$$

(ii) Determination of Na_2CO_3 .

1L of 1N HCl = 53g of NaOH

$$5 \text{ ml of } 0.101 \text{ N HCl} = 53 \times \frac{5 \times 0.101}{1000} \Rightarrow y_2 \text{ gram of Na}_2\text{CO}_3$$
$$= 0.027 \text{ g.}$$

$$\text{Strength of Na}_2\text{CO}_3 = y_2 \times 100 = 2.7 \text{ g/L}$$

Result \Rightarrow The given alkali mixture contains $\text{NaOH} = 4.6 \text{ g/L}$

The given alkali mixture contains $\text{Na}_2\text{CO}_3 = 2.7 \text{ g/L}$

Xurvi
Soni

Expt. No. _____

RESULTS - The given alkali mixture contains NaOH
 $= 4.6 \text{ g/L}$

The given alkali mix. contains $\text{Na}_2\text{CO}_3 = 2.7 \text{ g/L}$

Expected Clos / Daily life Applications

The total basic content of an antacid tablet can be determined in a similar manner.

PRECAUTION

- 1.) Rinse the pipette and burette.
- 2.) Always put the solution from burette drop wise to the conical flask.

Teacher's Signature : _____