
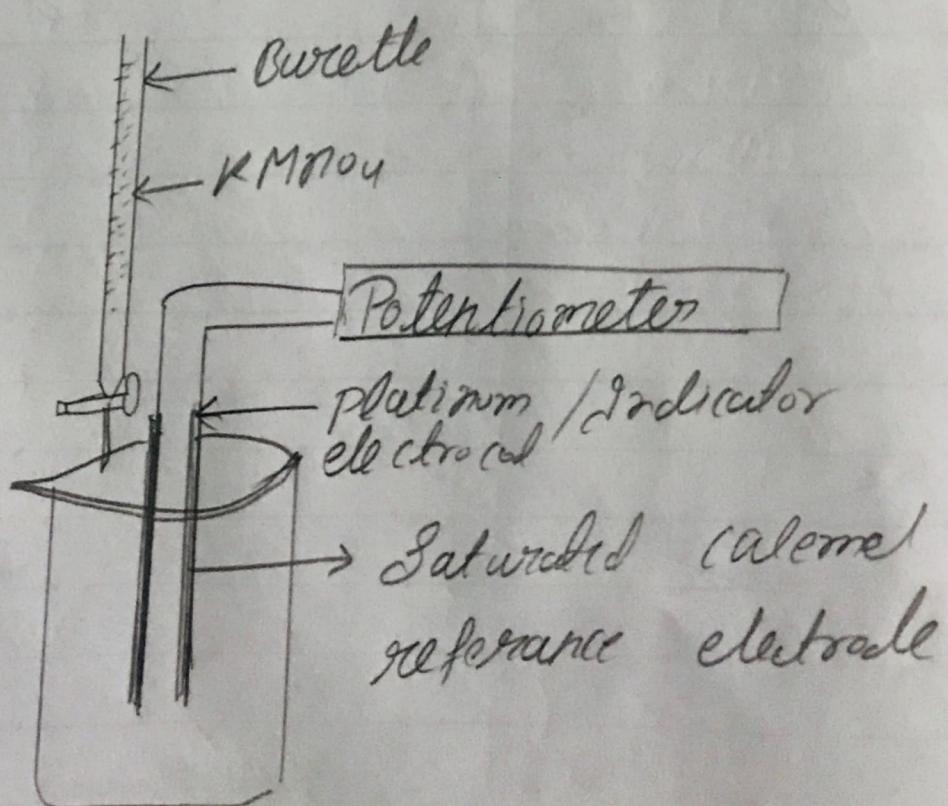
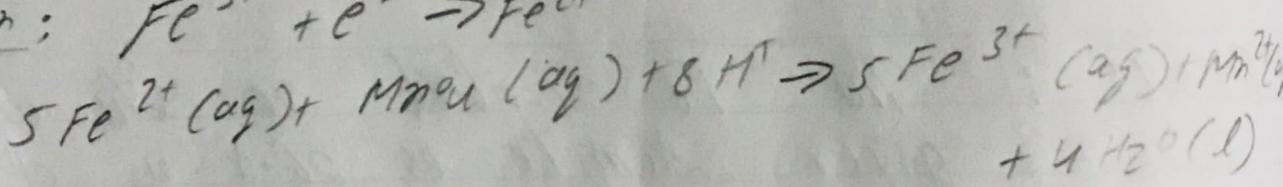


Experiment - 10

Experiment : To titrate potassium ferrocyanide against ammonium sulphate solution and determine the standard electrode potential of ferrous-ferric sulphate.

Apparatus : pipette, burette, beaker, funnel
(stand, clamp, potentiometer, calomel electrode, Ag / Hg electrode and platinum electrode).

Chemical : Mahr's salt solution, KMnO_4 and H_2SO_4

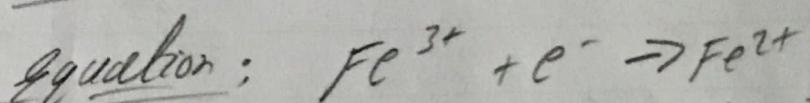

Experiment - 10

Experiment : To titrate potentiometrically ferrous ammonium sulphate solution against potassium permanganate and to determine the standard electrode potential of ferrous - ferroc system.

Apparatus : Pipette, burette, beaker, funnel, stand, clamp, galvanometer, calomel electrode (or, agc electrode) and platinum electrode.

Chemicals :- Mohr's salt solution (ferrous ammonium sulphate - $\text{FeSO}_4 \cdot (\text{NH}_4)_2\text{SO}_4 \cdot 6\text{H}_2\text{O}$) potassium permanganate (KMnO_4) and Sulphuric acid (H_2SO_4)


Theory : An electrochemical cell is a device which establishes measurable electrical potential differences and in which flow of electrical current is accompanied by an overall chemical reaction. A reversible cell is that in which the overall chemical reaction can be reversed on the presence of an opposing external electromotive force of magnitude greater than that of cell itself. An electrochemical

Observations :


(ii) Volume of $KMnO_4$ vs e.m.f of the solution
In beaker : FAS Solution ($0.01N$) + H_2SO_4 (IN)

$S \times 10$	Volume of $KMnO_4$ added (ml) (DDB applied)	EMF (volts)
1	0	0.300
2	0.2	0.341
3	0.7	0.361
4	1.2	0.372
5	1.7	0.380
6	2.2	0.388
7	2.7	0.394
8	3.2	0.400
9	3.7	0.405
10	4.2	0.410
11	4.7	0.415
12	5.2	0.419
13	5.7	0.423
14	6.2	0.427
15	6.7	0.432
16	7.2	0.436
17	7.7	0.441
18	8.2	0.446
19	8.7	0.451
20	9.2	0.456
21	9.7	0.464

cell consists of two electrodes on half cells, whose electrolytic solutions are either directly in contact with each other or connected through an intervening electrolytic solution. The net chemical change takes place at the individual electrodes, one of which is oxidation and the other reduction. At the reversible electrodes, the oxidized and reduced status of a system exist in equilibrium in the solution, where an inert metal electrode (like it) is dipped into it e.g., $\text{Fe}^{3+}/\text{Fe}^{2+}$ in which the reaction is:

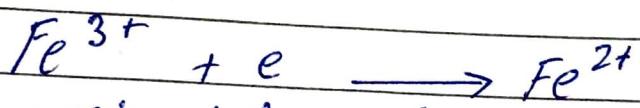


The experimental cell:

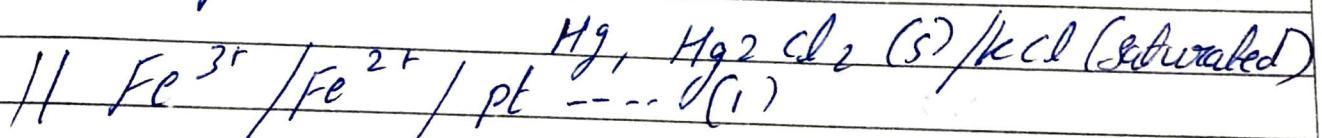

The emf of the cell:

$$E - E^\circ (\text{Fe}^{2+}/\text{Fe}^{3+}) + \frac{2.303RT}{nF} \log \frac{[\text{Fe}^{3+}]}{[\text{Fe}^{2+}]} - E$$

chemical reaction during polarimeter titration.

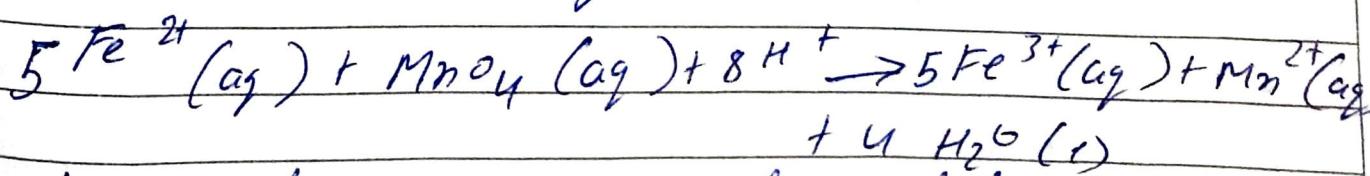

When potassium permanganate solution is added to Mohr's salt solution, the coagulation

S no	Volume of KMnO ₄ added ml	EMF (volts)
22	10.2	0.475
23	10.7	0.483
24	11.2	0.494
25	11.7	0.546
26	12.2	1.004
27	12.7	1.05
28	13.2	1.063
29	13.7	1.065
30	14.7	1.067
31	15.7	1.071
	16.7	

General calculation:

Emf at half equivalence point

$$\text{as } (A) = 0.423$$

$$E(Fe^{3+}/Fe^{2+}) = E(\text{half equivalent}) + E(\text{calm})$$

$$E(Fe^{3+}/Fe^{2+}) = A + 0.842 = 0.423 + 0.202 \\ = 0.665 \checkmark$$

$$E(Fe^{3+}/Fe^{2+}) = 0.665 \checkmark$$

of Fe^{3+} ions decreases and that the Fe^{2+} ions increases, and a result the unit of the cell increase slowly. Near to the equivalence point, an inflection is seen due to fall in concentration of Fe^{2+} ions ultimately to 0, resulting in sudden rise in current of the cell.

At

half equivalence, i.e. 2 drops

$$E(\text{half-equivalence}) = E^\circ(\text{Fe}^{3+}/\text{Fe}^{2+}) - E(\text{calomel})$$

$$= E^\circ[\text{Fe}^{3+}/\text{Fe}^{2+}] - 0.242 \text{ (at } 298\text{K)}$$


$$\therefore E^\circ(\text{Fe}^{3+}/\text{Fe}^{2+}) = E(\text{half equivalence}) + 0.242 \quad \dots \dots (3)$$

Procedure :- Take 50 ml of 0.1N FAS in beaker and add 5ml of 1N H_2SO_4 . Dip the platinum and standard calomel electrode in the solution.

2. Connect the indicator and reference electrodes to Bech and red terminals of potentiometer respectively.

3. Rinse and fill the burette with KMnO_4 (0.2N)

Teacher's Signature Ayush

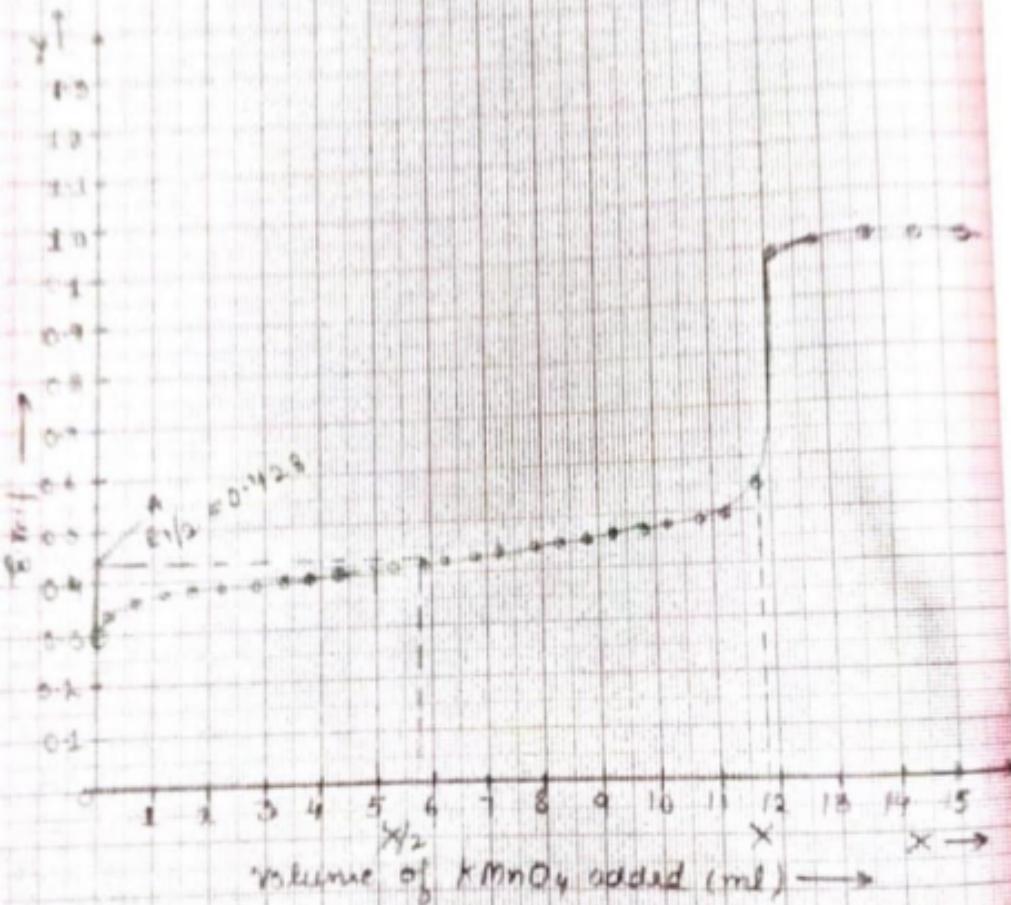

4 Note the initial inf of cell and start addition of titrant (KMnO₄) in portions of 1ml each. Near the equivalence point, pt. decrease the volume of additional titrant to 0.5 ml and later on to 0.2 ml and note the reading after each titration.

5 Continue to take 10-12 readings more, after a change in emf is notice.

6 Plot emf in volts, against the volume of KMnO₄ solution added (ml) and not sum the equivalence and potential at half the equivalence pt

General Calculation:

The emf at equivalence pt (A) observed from graph.

$$E(\text{half equivalence}) = E^{\circ} [\text{Fe}^{3+}/\text{Fe}^{2+}] - E(\text{calomel})$$

$$E[\text{Fe}^{3+}/\text{Fe}^{2+}] = E(\text{half quinacrine}) + E(\text{calomel})$$

Result

The standard half-cell potential of Fe²⁺/Fe³⁺ couple is 0.665 V

Teacher's Signature _____ Ayslfh.