
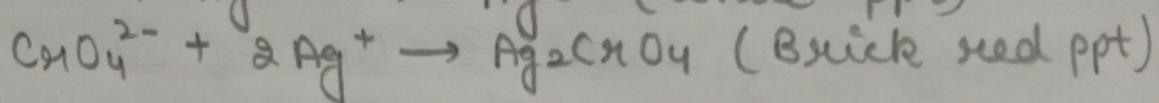
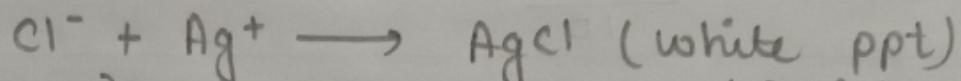
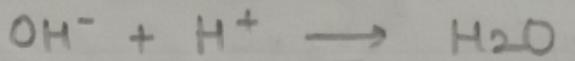
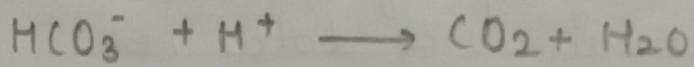


Experiment-6


Experiment: Find out the total alkalinity and chloride content in a water sample

Apparatus: Pipette, burette, beakers, conical flask, burette stand and clamp

Chemicals: Water sample, potassium chromate (K_2CrO_4), silver nitrate ($AgNO_3$), methyl orange and sulphuric acid (H_2SO_4),

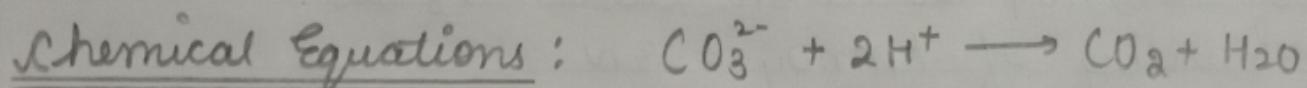

Structures:

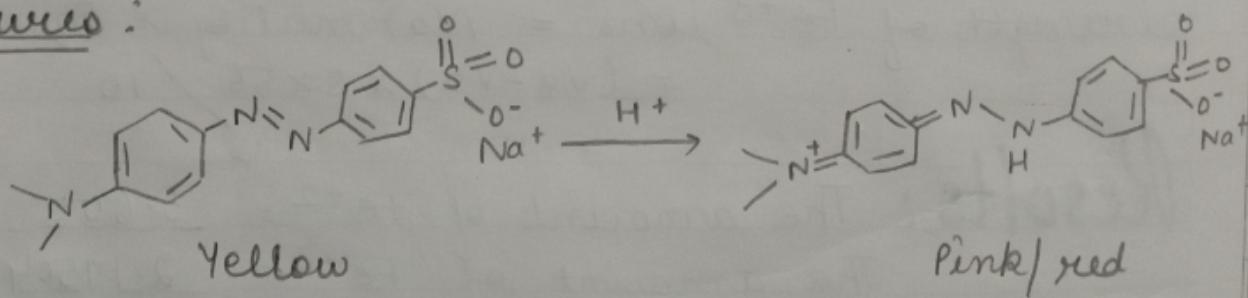

Indicator: Methyl orange

End point: Yellow to light pink

Procedure:

P.T.O.

EXPERIMENT-6

Experiment : Find out the total alkalinity and chloride content in a water sample

Apparatus : Pipette, burette, beakers, conical flask, burette, stand and clamp.

Chemicals : Water samples, potassium chromate (K_2CrO_4), silver nitrate ($AgNO_3$), methyl orange and sulphuric acid (H_2SO_4)

Theory : Alkalinity of water is due to the presence of hydroxides, carbonates and bicarbonates of the salts of calcium, magnesium, sodium and potassium. Similarly, the chloride content of water is due to the presence of chloride ions of these cations. Total alkalinity is estimated by titrating a known volume of water against a standard acid ($N/20 H_2SO_4$) using methyl orange as indicator in the neutral medium.



Chloride content is estimated by titrating a

Teacher's Signature

(i) Determination of total alkalinity of tap water The volume of tap water taken = 50 ml

Sr.no	Burette reading (ml)		Vol. of N/20 H ₂ SO ₄ used (ml)
	Initial	final	
1	0	6.9	6.9
2.	0	6.9	6.9
3.	0	6.9	6.9

Mean volume of H₂SO₄ used (V_1) = 6.9 ml

(ii) Determination of chloride contents of water sample Water sample in a titration flask = 10 ml

Sr.no	Burette reading (ml)		Vol. of N/100 AgNO ₃
	Initial	final	
1	0	12.9	12.9
2.	0	12.9	12.9
3.	0	12.9	12.9

Mean volume of AgNO₃ used V_2 = 12.9 ml

Calculations :


(i) Alkalinity

Applying the normality equation
(Tap water)

$$N_1 V_1 = N_2 V_2 \quad (\text{H}_2\text{SO}_4)$$

50 ml of tap water (of normality N_1) = x ml of $N/20 \text{ H}_2\text{SO}_4$

known volume against a standard silver nitrate solution ($N/100$) using potassium chromate as an indicator in the neutral medium

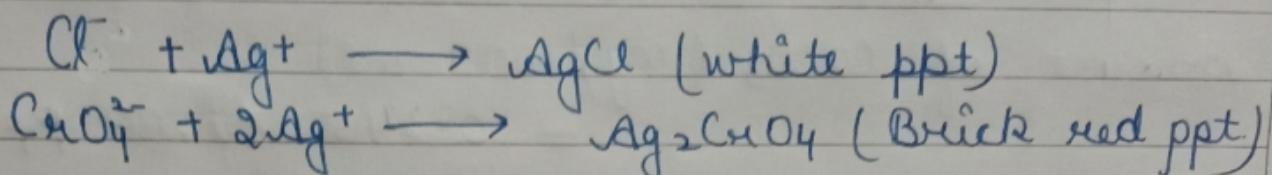

If water is found to be acidic, it is made neutral by adding solid calcium carbonate. In this case, some calcium carbonate must remain settled at the bottom. The results are expressed in parts per million (ppm)

Procedure :

(i) Determination of total alkalinity of tap water.

1. Wash, rinse and fill the burette with $N/20 \text{ H}_2\text{SO}_4$
2. Transfer .50 ml of tap water in the titration flask. Add 2-3 drops of methyl orange and titrate it against $N/20 \text{ H}_2\text{SO}_4$ till the colour changes from yellow to light pink, as an end point
3. Note the volume of the solution used and repeat the titration at least 5 times and

Teacher's Signature _____

$$V_2 = x = 6.9 \text{ ml}$$

Anjali Rana

$$N_1 = \frac{N_2 V_2}{V_1} = \frac{0.05x}{50} = \frac{0.05 \times 6.9}{50} = 0.0069$$

$$\begin{aligned} \text{Eq. wt of } \text{CaCO}_3 &= 50; & \text{Amount of } \text{CaCO}_3 (\text{gm/L}) \\ &= \text{Normality} \times \text{Eq. wt} \\ &= 0.0069 \times 50 = 0.345 \end{aligned}$$

$$\begin{aligned} \text{Amount of } \text{CaCO}_3 (\text{mg/1000 ml}) \\ &= 0.345 \times 1000 \text{ ppm} \end{aligned}$$

$$= 345 \text{ ppm}$$

(ii) Chloride content

Applying the normality equation,
(Tap water) $N_1 V_1 =$ (AgNO₃) $N_2 V_2$

$$10 \text{ ml of tap water (of normality } N_1) = y \text{ ml of } \frac{N/100}{\text{AgNO}_3 \text{ solution}}$$
$$y = 12.9 \text{ ml}$$

$$N_1 = \frac{N_2 V_2}{V_1} = \frac{y}{1000} = \frac{0.0129}{1000}$$

$$\text{Eq. wt of Cl}^- = 35.5$$

$$\begin{aligned} \text{Chloride content (gm/L)} &= \text{Normality} \times \text{Eq. wt} \\ &= \frac{y}{1000} \times 35.5 \\ &= 0.45795 \end{aligned}$$

$$\begin{aligned} \text{Chloride content (mg/1000 ml)} &= \frac{y \times 35.5 \text{ ppm}}{447.3 \text{ ppm}} \\ &= \frac{y}{12.3} \end{aligned}$$

take the mean of the closely related reading (x ml.)

(ii) Determination of chloride of water sample

1. Take 10 ml of water sample in a titration flask.
2. Add 3-4 drops of K_2CrO_4 and titrate against N/100 $AgNO_3$ from the burette till the appearance of light brick red color.
3. Note the volume of the used and repeat the titration at least 5 times and take the mean of the closely readings (y ml).

Observations.

let the volume of N/20 H_2SO_4 used = x ml

let the volume of N/100 $AgNO_3$ used = y ml

General Calculations

(i) Alkalinity

Applying the normality equation
(H_2SO_4)

$$N_1 V_1 = N_2 V_2$$

100 ml of tap water (of normality N_1) = x ml of
N/20 H_2SO_4

$$N_1 = \frac{N_2 V_2}{V_1} = \frac{0.05x}{100} = \frac{x}{2000}$$

Teacher's Signature _____

Results

Amount of total alkalinity in water sample 345 ppm of CaCO_3


Amount of chloride content in water sample 447.3 ppm

$$\text{Eq. wt. of } \text{CaCO}_3 = 50; \quad \begin{aligned} \text{Amount of } \text{CaCO}_3 (\text{gm/L}) \\ = \text{Normality} \times \text{Eq. wt.} \\ = \frac{x}{2000} \times 50 = \frac{x}{40} \end{aligned}$$

$$\begin{aligned} \text{Amount of } \text{CaCO}_3 (\text{mg/1000 mL}) \\ = \frac{x}{40} \times 1000 \text{ ppm} \end{aligned}$$

(ii) Chloride content:

Applying the normality equations

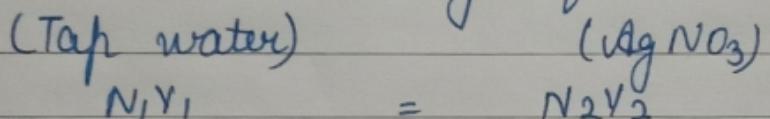

10 mL of tap water (of normality N_1) = y mL of $N/100$ AgNO_3 solution

$$N_1 = \frac{N_2 V_2}{V_1} = \frac{y}{1000}$$

Eq. wt. of Cl^- = 35.5

$$\begin{aligned} \text{Chloride content (gm/L)} &= \text{Normality} \times \text{Eq. wt.} \\ &= \frac{y}{1000} \times 35.5 \end{aligned}$$

$$\text{Chloride content (mg/1000 mL)} = y \times 35.5 \text{ ppm}$$

Results

Amount of total alkalinity in water sample 345 ppm of CaCO_3

Amount of chloride content in a sample 447.3 ppm